

#### 4. Preparation of TMB Substrate

Allow to equilibrate to room temperature. Ready for use.

#### 5. Reconstitution of Lyophilized Standards

Reconstitute one vial of each standard (Level I, II, III) by adding 0.15 ml of glass distilled or deionized water to each vial. Gently resuspend the pellet. Allow the rehydrated standards to stand at room temperature (15 to 30°C) for 10 minutes prior to use. **Do not** vortex or mix vigorously. Vigorous mixing will cause denaturation of the serum CIC and result in abnormally high values. **Reconstituted standards are stable for two working days. Store at 2 to 8°C when not in use.** Each reconstituted vial of standard is sufficient for 3 wells.  
NOTE: Reconstituted standards which have been diluted 1: 5 with Specimen diluent must be discarded after one working day.

#### 6. Reconstitution of Lyophilized Controls

Reconstitute one vial of each control by adding 0.25 ml of glass distilled or deionized water to each vial. Gently resuspend the pellet. Allow the rehydrated controls to stand at room temperature (15 to 30°C) for 10 minutes prior to use. **Do not** vortex or mix vigorously. Vigorous mixing will cause denaturation of the serum CIC and result in abnormally high values. **Reconstituted controls are stable for two working days. Store at 2 to 8°C when not in use.** Each reconstituted vial of control is sufficient for 6 wells.  
NOTE: Reconstituted controls which have been diluted 1: 5 with Specimen diluent must be discarded after one working day.

#### Test Procedure

- Using disposable dilution tubes, prepare 1:5 dilutions of the three reconstituted standards and the serum specimens in Specimen diluent. Mix gently. **Do not** vortex. Example: 0.05 ml of standard or serum specimen plus 0.20 ml of Specimen diluent (0.250 ml) is sufficient volume for duplicate testing.
- Take the microplate from the foil pouch by cutting above the ziplock. Identify the number of microplate wells required for the assay. Remove the strips (wells) from the holder not needed for the assay and replace into the pouch and reseal with the ziplock provided, leaving the desiccant packs inside. When stored at 2 to 8°C, resealed pouches will retain their reactivity. A pouch which has not been resealed or refrigerated after opening should be discarded.
- Wash the microwell strips, including Reagent Blank well (RB), three times with 1X Wash Buffer. After the last wash, tap microwell strips firmly on absorbent paper to remove excess Wash Buffer. Use the microwell strips immediately after washing. **Do not allow wells in the microwell strips to dry.**
- The upper left well (RB) on each frame is reserved for blanking the microwell strip reader. Refer to the manufacturer's instructions for blanking your reader.
- Mark the microwell strip holder to indicate location of Reagent Blank (RB), standards and serum specimens. NOTE: The standards and the serum specimens in steps

7-9 must be added to the microwell strips **within 10 minutes**. Pipette tips should be changed for each different standard and test specimen. Add 100 µl of each diluted standard and controls, in duplicate, to the appropriate wells.

- Add 100 µl of Specimen diluent to the reagent blank (RB) well.
- Add 100 µl of each diluted standard and control, in duplicate, to the appropriate wells.
- Add 100 µl of each diluted serum specimen, in duplicate, to the appropriate wells.
- Cover the microwell strip holder with a cover and incubate in a moist environment at 32 to 37°C for 60 minutes.
  - In a 37 C water bath: The covered microwell strips should not be floated on the water, but rested on a level support in the water bath, close to the surface of the water. Keep bath covered during incubation. Air temperature at level of microwell strips should be 32 to 37°C.
  - In a 37°C air incubator: Place a container with a cover inside the incubator. Leave this container inside the incubator to be sure the temperature inside is 37°C. Maintain the moisture by placing damp paper towels in the bottom of the container. Place microwell strip holder inside the container and keep the cover on for the incubation period.
- Prepare conjugate. Refer to Preparation Reagents.
- Wash the microwell strips, including blanking well(s), three times with 1X Wash Buffer. The washing procedure must be thorough. If an ELISA multichannel wash bottle is used:
  - Empty wells by vigorously shaking out contents.
  - Invert wash bottle and thoroughly flush wells with Wash Buffer by squeezing the bottle.
  - Remove Wash Buffer by vigorously shaking out contents.
  - Repeat steps 11b and 11c twice.
  - After the final wash, tap the inverted holder firmly on a clean paper towel to remove excess Wash Buffer.
  - Proceed immediately to the next step.

Gloves should be worn during washing steps. Household bleach should be added to waste material prior to disposal. NOTE: **Do not allow wells to dry during the test.** Drying of the wells may result in falsely high absorbance values.

- Using a multichannel micropipette, add 100 µl of conjugate to all wells, **including reagent blank well (RB)**.
- Incubate the microwell strips (uncovered) at room temperature in the dark for 60 minutes.
- Wash the microwell strips as described in Step 11. Proceed immediately to the next step.
- Using a multichannel micropipette, add 100 µl of substrate to all wells, **including the reagent blank well (RB)**.
- Incubate the microwell strips at room temperature in the dark for 30 minutes.

17. Using a multichannel micropipette, forcibly eject 100 µl of 1N sulfuric acid into each well, including the reagent blank well (RB). It is important that the sulfuric acid is spread quickly and uniformly throughout the well to completely inactivate the enzyme. Results can be read up to one hour after the sulfuric acid is added if the microwell strips are stored in the dark.

18. Set the microwell strip reader at a wavelength of 450 & 620 nm and blank the reader according to the manufacturer's instructions by using the reagent blank well (RB) containing substrate and acid.

19. Measure the color intensity in all microwell strips.

#### INTERPRETATION OF RESULTS

##### Standards and Controls

Calculate the mean absorbance value of the duplicate standards (Levels I, II, III) and controls. Duplicates should be within 20 percent of the mean.

##### Preparation of the Standard Curve and Determination of the Immune Complex Concentration (µg/ml HAG equivalents).

##### 1. Calculation Method


Calculations can be made on a variety of calculators supplied with linear regression programs. Follow the manufacturer's instructions for use of these programs to establish the standard curve and the immune complex concentration (µg/ml HAG equivalents) of each serum specimen.

In order for the assay to be valid, the following criteria must be met:

- Slope 0.008-0.013
- Y-intercept less than 0.250
- Correlation Coefficient (r<sup>2</sup>): greater than 0.950

##### 2. Graph Method

- Plot the mean absorbance values for standards (Level I, II, III) on the ordinate against the corresponding µg/ml HAG on the abscissa, using linear graph paper. Refer to vial labels of standards for assigned values expressed as µg/ml.
- Draw a straight line for the **best fit** of the plotted points.

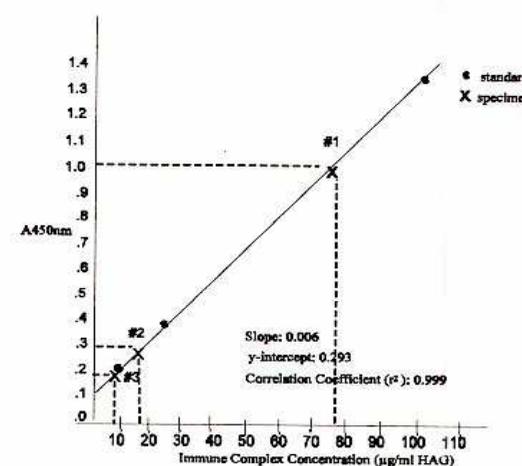

- Locate the point corresponding to the average absorbance value of each sera specimen and determine the corresponding µg/ml HAG equivalents.
- Record the immune complex concentrations for each serum as µg/ml HAG equivalents.

#### 3. Illustration of Standard Curve

A representative Standard Curve is shown in Figure 1. This curve cannot be used to derive test results. Each laboratory must prepare a Standard Curve for each group of microwell strips assayed (up to twelve). Refer to each vial label of Standard (Level I, II, III) for assigned values expressed as µg/ml.

Figure 1. Standard Curve

	Absorbance Value	Immune Complex Concentration (µg/ml HAG equivalents)
Standard Level I	0.155	6
Standard Level II	0.405	26
Standard Level III	1.368	106
Specimen 1	1.027	77
Specimen 2	0.294	17
Specimen 3	0.123	3.0

#### Criteria for Interpretation of Results

##### 1. Negative Results

Results from 0 to 8 µg/ml HAG equivalents can be considered negative for significant levels of CIC.

##### 2. Positive Results

Results of 9 µg/ml HAG equivalents or greater can be interpreted as positive for significant levels of CIC.

## LIMITATIONS OF PROCEDURES

- Falsely elevated results will be obtained if:**
  - Plasma is used instead of serum.
  - Specimens are collected as plasma in anticoagulants and subsequently defibrinated.
  - Serum specimens have been heat inactivated.
  - Specimens are improperly handled or stored.
- Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results.
- Adherence to the humidity, temperature and time periods for incubation is essential for accurate results.
- Bubbles or excess reagent in the wells may result in cross contamination of the wells when a microwell cover is used.
- Bacterial or fungal contamination of either serum specimens or reagents or cross contamination between reagents may cause erroneous results.
- A small percentage of the population will have values close to the statistically derived upper limit of normal. These results should be interpreted as borderline. For these individuals, retesting is recommended.
- If another manufacturer's test is used to determine CIC, the laboratory should establish normal ranges for the test being used.
- CIC formation may be associated with pathological conditions other than RA and SLE.


## SPECIFIC PERFORMANCE CHARACTERISTICS

### Accuracy

A World Health Organization (WHO) standard for immune complex assays containing 800 micrograms of heat aggregated IgG per ml was used to prepare standards in fresh frozen human serum ranging in concentrations from 10 to 150 µg/ml. The samples were coded and assayed as unknowns using **C3d SCANLISA** Test System for Circulating Immune Complexes (C3d-IgG). The assayed values are shown plotted on the abscissa against the expected values plotted on the ordinate in Figure 2. The actual assay values appear in parentheses. Each point represents the mean of triplicate absorbance values. The response was assumed linear over the entire test range and gave a coefficient of determination of 0.995.

Assay of WHO Standard for Immune Complex Assays using **C3d SCANLISA** Test System for Circulating Immune Complexes (C3d-IgG).

Figure 2.

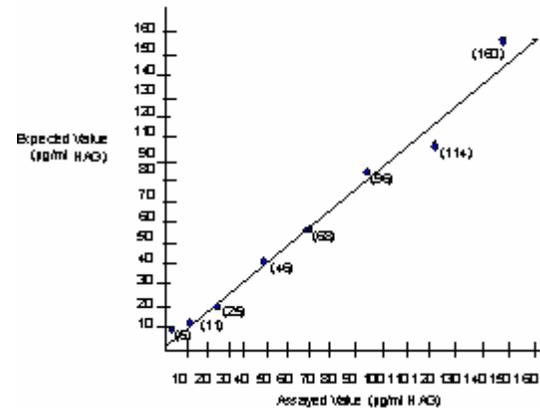

### Precision

Six serum specimens were tested in duplicate on two different microwell strips on the same day and converted to µg/ml HAG equivalents. The coefficient of variation between strips was less than 5%. Two hundred sixty-four replicates of a single serum specimen were assayed on each of twenty-two microwell strips. The average coefficient of variation of the four plates was less than 5%.

### BIBLIOGRAPHY

- Theofilopoulos AN, Dixon FK. The biology and detection of immune complexes. In: Dixon FJ, Kunkle HG; eds. *Advances in immunology*, New York; Academic Press; 1979; 28: 89-220.
- Nydegger UE, Davis JS, IV. Soluble immune complexes in human disease, *CRC Crit Rev Clin Lab Sci* 1980; 12(2): 123-170.
- Williams RC, Jr. Immune complexes in human diseases. *Ann Rev Med* 1981; 32: 13-28.
- Nydegger UE, Lambert PH, Gerber H, Miescher PA. Circulating immune complexes in the serum in systemic lupus erythematosus and in carriers of hepatitis B antigen. *J Clin Invest* 1974; 54: 297-309.
- Lambert, PH, et al. Immune complexes and complement activation in rheumatoid arthritis. *Proceedings of the symposium on immunology and infection in the rheumatic diseases*. London: Blackwell and Co. Ltd. 1974.
- Theofilopoulos AN, Wilson CB, Dixon FK. The Raji cell radioimmune assay for detecting immune complexes in human sera. *J. Clin Invest* 1976; 57:169.
- Pereira AB, Theofilopoulos AN, Dixon FJ. Detection and partial characterization of circulating immune complexes with solid-phase anti-C3. *J. Immunol* 1980; 125:763-770.

5. Do **not** use plasma. The serum should be separated from the clot within eight hours of collection and transferred to a plastic test tube. The test specimen can be centrifuged at 850-1500 xg for 10 minutes at ambient temperature (15 to 30°C). Mechanical serum separated devices such as ISO Filter† or equivalent may be used. The resulting test specimen must be completely free of cellular debris. Once the serum has been separated from the clot, it may be stored at 2 to 8°C for up to 24 hours.

If the serum specimen cannot be tested within one day of collection, the specimen must be stored at -70°C until tested. Serum specimens are stable at -70°C for at least three months. Freeze and thaw only once. Specimens may be thawed either at room temperature or for up to 5 minutes at 37° C.

### Materials Required But Not Provided

- 1 ml, 5 ml and 10 ml serological pipettes
- 100 µl adjustable single-channel micropipette
- 100 µL multichannel micropipette
- 50 µl to 200 µl disposable pipette tips (recommended by multichannel micropipette manufacturer)
- Multichannel micropipette reservoir
- Disposable 12x75 mm or 13x100 mm glass tubes
- Disposable plastic 25 ml tube with cap
- Plastic microwell strip cover or polyethylene film
- ELISA multichannel wash bottle or automatic (semi-automatic) washing system
- Aliquot mixer
- 37°C humidified incubator or water bath with cover
- Microwell strip reader capable of reading at 450 nm
- Glass-distilled or deionized water

### Precautions

- Disposable glass or plastic must be used for preparation of the conjugate.
- Cross contamination between reagents may invalidate the test results. The reagents should be visually inspected for evidence of bacterial or fungal contamination. Permanently labeled, dedicated multichannel micropipette reservoirs for the appropriate reagents are recommended.
- Handle the TMB with care. Avoid contact with eyes, skin or clothing as TMB may cause irritation or an allergic skin reaction.
- The TMB is light sensitive. Keep vial **tightly** closed when not in use.
- All micropipettes must be used with care and calibrated regularly following the manufacturer's instructions.
- The microwell strips must be thoroughly washed during the washing steps.
- All serum specimens should be handled as if capable of transmitting disease.
- Wash Buffer, Specimen diluent and distilled or deionized water must be a room temperature prior to use.

- Diluted serum specimens and standards should be added to the microwell strips at room temperature (15 to 30°C).
- During the second and third (last) washing procedure, disposable gloves should be worn and all reagent and washing solutions should be discarded into an appropriate container. It is recommended that household bleach be added to the waste material prior to disposal. **Do not** add bleach prior to washing because of possible splash back.

## Preparation of Reagents

- Wash Buffer (10X):** Dilute the buffer 1 + 9 using distilled or deionized water or by pouring total contents into a graduated cylinder and add distilled or deionized water to 1,000 ml final volume. Mix thoroughly and store the 1X solution at 2 - 8°C.
- Specimen diluent** is "Ready to Use". DO NOT Dilute further. Allow to equilibrate to room temperature.
- Preparation of Conjugate**  
Disposable glass or plastic must be used. The Conjugate Concentrate must be diluted with Specimen diluent just prior to use. The dilution ratio for the Conjugate Concentrate appears on the label. For example, if the label states dilution ratio 1:100: then add 0.12 ml of Conjugate Concentrate to 11.86 ml Specimen diluent. A volume of 12 ml is needed for twelve (12) microwell strips.

The conjugate should be prepared just prior to use and should be used within 30 minutes.

Number of Strips	Conjugate			
	Ratio	Total mL.	Conj. Conc. (µl)	Specimen diluent. (ml)
1-3	1:100	3.0	30	2.97
	1:120	3.6	30	3.57
	1:140	4.2	30	4.17
4-6	1:100	6.0	60	5.94
	1:120	7.2	60	7.14
	1:140	8.4	60	8.34
7-9	1:100	9.0	90	8.91
	1:120	10.8	90	10.71
	1:140	12.6	90	12.51
10-12	1:100	12.0	120	11.88
	1:120	14.4	120	14.28
	1:140	16.8	120	16.68
Stability at RT:		30 minutes		

† Trademark Becton Dickenson

## SUMMARY AND EXPLANATION

The harmful effects of circulating immune complexes (CIC) have been thoroughly reviewed by Theofilopoulos and Dixon and others.<sup>1-3</sup> These CIC are not normally expressed in the serum of healthy individuals but are readily detected in the serum of patients with rheumatoid arthritis (RA) and systemic lupus erythematosus (SLE).<sup>4,5</sup> CIC are particularly prominent in RA and SLE patients during active disease.

No single procedure appears to detect all types of CIC, but those procedures such as the Raji Assay and C3d Solid Phase Assay which detect CIC containing fragments of C3 (i.e., C3d) appear to detect clinically relevant events.

C3d SCANLISA Strip, ELISA Test System for Circulating Immune Complexes (C3d-IgG) detects CIC containing both C3d and IgG.<sup>6</sup> A positive result is indicated by the presence of CIC levels equal to or greater than 9 µg/ml. A negative result is indicated by a level CIC less than 9 µg/ml. The concentration of CIC is measured as µg/ml heat aggregated human globulin (HAG) equivalents.

## PRINCIPAL OF THE PROCEDURE

The procedure is a three stage test carried out on microwell strips which have been coated with monoclonal antibody specific for the C3d component of complement (solid phase).

In the first stage, a diluted test serum is added to a test well and incubated for a specified length of time. If CIC containing C3d are present in the test serum, they will attach to the solid phase anti-C3d. If CIC are not present, all serum components will be removed in the subsequent washing step.

In the second stage, monoclonal anti-human IgG conjugated to horseradish peroxidase enzyme (conjugate) is added to the test well. The conjugate will bind specifically to the IgG of the bound CIC. If CIC are not present in the serum, the unbound conjugate will be removed in the subsequent washing step.

In the third stage, enzyme substrate is added to the test well. If bound conjugate is present the substrate will be oxidized, resulting in a colored end product. Acid is added to stop the reaction and fix the color. The color intensity is proportional to the amount of bound conjugate and therefore to the amount of CIC present in the serum specimen. The color intensity is measured with a microwell reader.

## REAGENTS

Reagents are standardized by lot number and must not be interchanged with reagents from a kit with a different lot number.

**Kit Components (96 Tests),** Catalog # C3196

1 Plate (12 X 8 break-apart wells) coated with Monoclonal Antibody (Murine) to Human C3d in a ziplock foil pouch.

1 Bottle 10X Wash Buffer Concentrate – dilute 1:10

2 Bottles Specimen diluent – “Ready to Use”

1 Vial Conjugate Concentrate—Monoclonal anti-human IgG (Murine) conjugated to horseradish peroxidase; contains 0.1% Proclin 300®

3 Vials Standard Level I (lyophilized human serum)

3 Vials Standard Level II (lyophilized human serum)

3 Vials Standard Level III (lyophilized human serum)

2 Vials Positive Control (lyophilized human serum)

2 Vials Negative Control (lyophilized human serum)

1 Vial Tetramethylbenzidine (TMB) Substrate “Ready to Use”

1 Bottle Stop Solution (1N H<sub>2</sub>SO<sub>4</sub>) – “Ready to Use”

Instructions for preparation of reagents are provided in the section entitled: PROCEDURE—Preparation of Reagents.

For in vitro diagnostic use.

Store at 2 to 8°C. Do not freeze.

Source material from which reagents of human origin were derived was found nonreactive for HBsAg and HTLV-III when tested with licensed reagents. No known test method can offer assurance that products derived from human blood will not transmit hepatitis or other infectious agents.

**CAUTION: HANDLE AS IF CAPABLE OF TRANSMITTING INFECTIOUS AGENTS.**

## SPECIMEN COLLECTION

No special preparation of the individual is required prior to specimen collection. The whole-blood test specimen should be collected by accepted medical techniques. **The specimen must be allowed to clot for two hours at room temperature (15 to 30°C) before the serum is removed.**

1. Do **not** allow the specimen to clot at refrigerator temperature (2 to 8°C).
2. Do **not** store the clotted specimen at 2 to 8°C.
3. Do **not** use anticoagulated specimens collected as plasma and then clotted.
4. Do **not** heat inactivate serum specimens.

**C3d SCANLISA  
Test System for  
Circulating Immune Complexes  
(C3d-IgG)**

**FOR *IN VITRO* DIAGNOSTIC USE**

CATALOG NO.:

C3196

96 Test

**New Phase II**

***\* Ready to Use Specimen diluent,  
Stop Solution, 10X Wash Buffer,  
and 100  $\mu$ L / well Sample Size***

 SCIMEDX CORPORATION  
400 Ford Road  
Denville, NJ 07834 USA  
Tel: 800.221.5598 Tel: 973.625.8822  
Fax: 973.625.8796  
www.scimedx.com

EC	REP	Emergo Europe Molenstraat 15 2513 BH The Hague The Netherlands
----	-----	---



Printed in U.S.A. Rev. C 12/09/03

C3196.C



*Your Preferred  
Autoimmune and  
Infectious Disease  
Detection Company*

400 FORD ROAD  
DENVILLE, NEW JERSEY 07834 U.S.A.  
(973) 625-8822 ♦ FAX: (973) 625-8796

CUSTOMER INQUIRIES:  
(800) 221-5598

